

**This Page Is Inserted by IFW Operations
and is not a part of the Official Record**

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images may include (but are not limited to):

- **BLACK BORDERS**
- **TEXT CUT OFF AT TOP, BOTTOM OR SIDES**
- **FADED TEXT**
- **ILLEGIBLE TEXT**
- **SKEWED/SLANTED IMAGES**
- **COLORED PHOTOS**
- **BLACK OR VERY BLACK AND WHITE DARK PHOTOS**
- **GRAY SCALE DOCUMENTS**

IMAGES ARE BEST AVAILABLE COPY.

**As rescanning documents *will not* correct images,
please do not report the images to the
Image Problem Mailbox.**

12

EUROPEAN PATENT APPLICATION

21 Application number: 86118191.5

5 Int. Cl. 4: C12Q 1/68

22 Date of filing: 01.12.86

43 Date of publication of application:
08.06.88 Bulletin 88/23

64 Designated Contracting States:
AT BE CH DE ES FR GB GR IT LI LU NL SE

71 Applicant: MOLECULAR BIOSYSTEMS, INC.
11180 Roselle Street, Suite A
San Diego, CA 92121(US)

72 Inventor: Heller, Michael J.
1357 Acton Avenue
Poway, CA 92064(US)

74 Representative: Patentanwälte Grünecker,
Kinkeldey, Stockmalr & Partner
Maximilianstrasse 58
D-8000 München 22(DE)

54 Method for increasing the sensitivity of nucleic acid hybridization assays.

57 An improvement in the sensitivity of hybridization assays for detecting polynucleotide target sequences is provided by selective dehybridization of the bound portion of the probe, separation of the dehybridized portion from the support and the residual sample, and concentrating the separated probe before detecting its presence by means of its reporter group. The concentration step can be carried out by adsorption of the dehybridized probe onto a basic ion exchange material. A displacer sequence having greater binding affinity for the target polynucleotide sequence may be used to dehybridize the bound probe portion for subsequent concentration.

EP 0 269 764 A1

METHOD FOR INCREASING THE SENSITIVITY OF NUCLEIC ACID HYBRIDIZATION ASSAYS

FIELD OF INVENTION, BACKGROUND AND PRIOR ART

The field of this invention involves a hybridization assay system for detecting target polynucleotide sequences in deoxyribonucleic acid (DNA) and/or ribonucleic acid (RNA) samples. More specifically, the invention is concerned with a method for increasing the sensitivity and where necessary the selectivity of the hybridization assay using labelled oligonucleotide probes to detect specific target sequences.

Hybridization assays may be employed for the detection and identification of DNA or RNA sequences. Published methods as used particularly in recombinant DNA research are described in Methods in Enzymology, Vol. 68, pp. 379-469 (1979); and Vol. 65, Part 1, pp. 468-478 (1968). One such method involving a preliminary separation of the nucleic acid fragments by electrophoresis is known as the "Southern Blot Filter Hybridization Method." See Southern, J. Mol. Biol. (1975) 98:503. Hybridization probe methods for detecting pathogens are described in U.S. Patent 4,358,539.

There has been a recognized need for improving hybridization assays. Instead of radiolabelled probes, which can require several days for development of the X-ray film, detection by means of luminescent-labelled probes are claimed to have advantages. (See published European Patent Application 0 070 687.) A procedure for obtaining increased sensitivity from luminescent-labelled probes is disclosed in published European Patent Application 0 070 685. Two probes are employed, having respectively a chemiluminescent catalyst and an absorber/emitter (fluorophore) moiety. Upon hybridization to the appropriate polynucleotide target sequence, the light signal from the probes is shifted to a higher wavelength. This wavelength shifting leads to an improved and more sensitive assay system.

As an alternative to employing long double-stranded probes resulting from recombinant DNA cloning procedures, short synthetically produced single-stranded probes have been described. (See, for example, PCT Publication Number WO 84/03285). Such short probes, which normally contain from 18 to 40 nucleotides, are available from Molecular Biosystems, Inc., San Diego, California. Typically, such short probes can be labelled with biotin, enzymes, and fluorescent reporter groups. These probes can be constructed complementary to a unique sequence of the target DNA or RNA fragments. Furthermore, they are in stable single-stranded form and renaturation cannot interfere with the assay. Short labelled DNA probes as described above, in principle, should be useful in developing rapid and convenient non-radioisotope DNA hybridization assays. With existing procedures, concentrations of the hybridized labelled probe can be detected down to about 10^{-16} to 10^{-17} moles under ideal conditions. For most of the potentially important clinical applications of DNA probes (infectious diseases, latent virus, genetic disorders, etc.), it will be necessary to have sensitivity levels of at least 10^{-18} moles and lower. With present techniques the inherent sensitivities possible with non-radioisotope labels such as fluorophores and enzymes are still limited by a wide variety of problems. In the case of fluorescent labels, background fluorescence and light scattering from the hybridization support materials (nitrocellulose and nylon filters, polystyrene beads, etc.) lead to significant loss of inherent sensitivity. In the case of enzyme labels, non-specific binding of the labeled probe and interfering stains and colors from sample material can limit ultimate sensitivity of this label. Other problems superimposed upon those mentioned above lead to further losses of sensitivity and limit the practicality of these systems for clinical diagnostic purposes.

SUMMARY OF INVENTION

The present invention provides a method for ultrasensitive detection of reporter group-labelled DNA probes used in hybridization assay systems. The method is particularly applicable to fluorescently-labelled short probes, but should also provide improvements for enzyme-labelled probes. In the case of fluorescent probes, greater sensitivities can be obtained than with any prior hybridization assay system, viz. as little as 10^{-18} moles of probe and below can be detected.

The method of this invention utilizes a probe concentration or refocusing procedure. After the basic hybridization procedure has been carried out and the support matrix (nitrocellulose, for example) has been washed free of unbound probe, instead of attempting to detect the presence of bound probe directly as in conventional systems, the probe is subjected to a dehybridization procedure. Dehybridization can be carried out by treatment with denaturing agents (heat, urea, etc.) which will release the bound probe from the support matrix and liberate it into the aqueous solution. The bound probe can also be liberated in a much

more selective fashion using a second DNA probe which is called a "displacer probe." In either case the liberated probe, in solution, is now concentrated or refocused by addition of a small amount of an appropriate material with high affinity for the probe.

In one preferred procedure, after dehybridization from the support matrix (nitrocellulose or nylon filters), the fluorescently labelled probe is adsorbed onto a small amount of adsorbent beads or particles having high affinity for oligonucleotides. For example, the adsorbent may comprise beads of porous, crosslinked dextran or glass particles substituted with positively charged groups. Such adsorbents are available commercially and are improved substrates for fluorescent probe detection in particular. For example, anion exchange materials such as diethylaminoethyl (DEAE) or quaternary aminoethyl (QAE) dextran and glass beads have extremely high affinity for labelled oligonucleotide probes. This high affinity is due to the electrostatic or ion exchange interaction between the polyanionic nucleotide and the high positive charge density on the beads. These beads have been found to provide extremely low background fluorescence and significantly less light scattering effects than nitrocellulose or nylon filter material. Most importantly, the dehybridization and adsorption onto a small amount of adsorbent material concentrates (10-5000 fold) the fluorescent probe onto a substrate much more favorable for fluorescent (epifluorescent) detection. In summary, the fluorescent signal is now highly concentrated, background is significantly reduced, and interferences from the hybridization support material (nitrocellulose filters, other sample components on the filters, etc.) have been eliminated.

While the specific procedures discussed above are preferred, within the broader context of the improved method of this invention other concentration or focusing procedures can be employed. Moreover, the probes can be used which are labelled with other types of reporter groups (i.e., enzymes or luminescent groups). With such modifications, improvements in detection sensitivity can also be obtained in accordance with the principles of this invention.

DETAILED DESCRIPTION

In practicing the method of this invention, the DNA or RNA test samples containing target sequences may be prepared by any one of a number of known procedures and attached to suitable immobilization support matrix. For example, such procedures are described in Methods in Enzymology, Vol. 68, pp. 379-469 (1979); and Vol. 65, Part 1, pp. 468-478 (1980), U.S. Patent 4,358,539; and in published European Patent Applications Nos. 0 070 685 and 0 070 687. Any of the usual supports can be employed, such as nitrocellulose filters, Zetabind® (nylon) filters, polystyrene beads, etc.

The following is a general description of the basic procedures for carrying out the hybridization of labelled short DNA probes to samples of target DNA immobilized on nitrocellulose and Zetabind® filters. With minor variations the procedure is applicable for hybridizing radioisotope (³²P), biotin, enzyme, lumiphore, and fluorophore labelled short probes.

Buffers and Solutions

All buffers should be made with double-distilled water. Concentrated stock solutions should be diluted with double-distilled water. Filter, where indicated, through 0.45μ nitrocellulose.

20 SSPE BUFFER. Dissolve 210 g NaCl, 27.6 g NaH₂PO₄ · H₂O, 7.4 g Na₂EDTA, and approximately 8.3 ml of saturated NaOH (to pH 7.0). Adjust volume to one liter, filter, autoclave and store at room temperature.

20 SSC BUFFER. Dissolve 175.3 g NaCl and 88.2 g sodium citrate in 800 ml water. Adjust pH to 7.0 with a few drops of concentrated NaOH. Adjust volume to one liter, filter, autoclave, and store at room temperature.

0.3 M NaOH. Store at room temperature.

1 M and 2 M Ammonium Acetate. Filter and store at room temperature.

4 x BP. Dissolve 2 g each bovine serum albumin, crystalline Pentex Fraction V., and polyvinylpyrrolidone, average mol. wt. = 40,000 in 100 ml water. Filter and store at -20°C in small aliquots.

10% (W/V) Sodium Dodecyl Sulfate (SDS). Store at room temperature.

Other Materials:

Nitrocellulose filters (Schleicher & Schuell, Keene, NH, #BA85) or Zetabind membranes (AMF, Meriden, CT, #NM511-01-04-5SP).

5 Whatman 3MM chromatography paper.

DEAE-cellulose. Whatman DE52, precycle according to the manufacturer.

Dot blot filtration apparatus (HYBRI-⁺ DOTTM, BRL, Bethesda, MD; Minifold II, S & S, Keene, NH; BIODOTTM, Bio-Rad Laboratories, Richmond, CA.)

Water baths.

10 Vacuum pump.

Vacuum oven.

Automatic micropipets and tips.

Disposable tuberculin syringes (1 ml).

Glass wool.

15

Sample Preparation

The method of fixing target DNA to nitrocellulose filters will vary depending on the source of the sample. Purified DNA may be bound by filtration by Southern blotting, by Dot blotting. Virus-infected whole cells can be fixed to nitrocellulose by spotting followed by denaturation.

The following procedure is suitable for binding the DNA to nitrocellulose filters. It should be modified according to the needs of the individual experiment. The following steps are performed at room temperature unless indicated otherwise.

25 NOTE: Nitrocellulose and nylon filters should be stored in a dry cool place and protected from fumes of organic solvents. Handle the filters with gloves.

1. Cut the nitrocellulose or nylon to fit the dot blot apparatus and wet the filters in hot (approximately 60°C) sterile water. Make certain that the sheet is completely wet, then rinse it several additional times with hot clean water.

30 2. Soak the sheet in 6 x SSC buffer while preparing the samples.

3. Resuspend DNA in 0.1 ml 0.3 M NaOH and boil for 5 minutes.

4. Immediately chill the sample in ice and add 0.1 ml 2 M Ammonium Acetate. Mix well.

5. Clamp sheet in the dot blot apparatus, apply a gentle vacuum and wash the wells with 200µl 6 x SSC.

35 6. Decrease the vacuum and add the sample. A sample should take several minutes to filter in order for the DNA to bind quantitatively.

7. After the samples have filtered through, wash the wells as in Step 5.

8. Remove the filter from the apparatus, wash it in 6 x SSC, and air dry.

40 9. Place the filter between the folded sheets of Whatman 3MM paper and bake in a vacuum oven for one hour at 80°C. Nylon support should be dried for 15-30 minutes at 80°C.

10. Filter sheets prepared thus may be stored dry at 4°C in sealed plastic bags.

Hybridization Detection

45

1. Prehybridize the filter in 5 x SSPE, 1 x BP, and 1% SDS for 5 minutes at 50° ± 2°C.

2. Transfer the filter to a boilable polyethylene bag, add 50µl per cm² of 5 x SSPE, 1 x BP, 1% SDS, containing ~ 5-500 ng/ml of appropriately labelled probe and heat-seal the bag. Incubate for 30 minutes at 50° ± 2°C with gentle agitation.

50 3. After the incubation, open the bag and discard the buffer into an appropriate waste container. Wash the filter as follows, using a large excess of 1 x SSPE, 1% SDS; three times at 37°C for five minutes each followed by one wash at 50°C for one minute.

4. At this point different procedures can be carried out depending upon the nature of the label and detection system (radio isotope, enzyme/color dye system, etc.). With regard to the further teaching in this patent application and in reference to fluorescently labeled probes in particular, the filters would be washed one time with 1 x SSPE. At this point procedures discussed later in the application apply.

The probe is designed and selected to provide a complementary sequence for binding to the target DNA or RNA sequences. Short-chain probes are preferred, such as those containing from 10 to 50

Other Materials:

Nitrocellulose filters (Schleicher & Schuell, Keene, NH, #BA85) or Zetabind membranes (AMF, Meriden, CT, #NM511-01-04-5SP).

5 Whatman 3MM chromatography paper.

DEAE-cellulose. Whatman DE52, precycle according to the manufacturer.

Dot blot filtration apparatus (HYBRIDOT™, BRL, Bethesda, MD; Minifold II, S & S, Keene, NH; BIODOT™, Bio-Rad Laboratories, Richmond, CA.)

Water baths.

10 Vacuum pump.

Vacuum oven.

Automatic micropipets and tips.

Disposable tuberculin syringes (1 ml).

Glass wool.

15

Sample Preparation

The method of fixing target DNA to nitrocellulose filters will vary depending on the source of the sample. Purified DNA may be bound by filtration by Southern blotting, by Dot blotting. Virus-infected whole cells can be fixed to nitrocellulose by spotting followed by denaturation.

The following procedure is suitable for binding the DNA to nitrocellulose filters. It should be modified according to the needs of the individual experiment. The following steps are performed at room temperature unless indicated otherwise.

25 NOTE: Nitrocellulose and nylon filters should be stored in a dry cool place and protected from fumes of organic solvents. Handle the filters with gloves.

1. Cut the nitrocellulose or nylon to fit the dot blot apparatus and wet the filters in hot (approximately 60°C) sterile water. Make certain that the sheet is completely wet, then rinse it several additional times with hot clean water.

30 2. Soak the sheet in 6 x SSC buffer while preparing the samples.

3. Resuspend DNA in 0.1 ml 0.3 M NaOH and boil for 5 minutes.

4. Immediately chill the sample in ice and add 0.1 ml 2 M Ammonium Acetate. Mix well.

5. Clamp sheet in the dot blot apparatus, apply a gentle vacuum and wash the wells with 200µl 6 x SSC.

35 6. Decrease the vacuum and add the sample. A sample should take several minutes to filter in order for the DNA to bind quantitatively.

7. After the samples have filtered through, wash the wells as in Step 5.

8. Remove the filter from the apparatus, wash it in 6 x SSC, and air dry.

40 9. Place the filter between the folded sheets of Whatman 3MM paper and bake in a vacuum oven for one hour at 80°C. Nylon support should be dried for 15-30 minutes at 80°C.

10. Filter sheets prepared thus may be stored dry at 4°C in sealed plastic bags.

Hybridization Detection

45

1. Prehybridize the filter in 5 x SSPE, 1 x BP, and 1% SDS for 5 minutes at 50° ± 2°C.

2. Transfer the filter to a boilable polyethylene bag, add 50µl per cm² of 5 x SSPE, 1 x BP, 1% SDS, containing ~ 5-500 ng/ml of appropriately labelled probe and heat-seal the bag. Incubate for 30 minutes at 50° ± 2°C with gentle agitation.

50 3. After the incubation, open the bag and discard the buffer into an appropriate waste container. Wash the filter as follows, using a large excess of 1 x SSPE, 1% SDS; three times at 37°C for five minutes each followed by one wash at 50°C for one minute.

4. At this point different procedures can be carried out depending upon the nature of the label and detection system (radio isotope, enzyme/color dye system, etc.). With regard to the further teaching in this patent application and in reference to fluorescently labeled probes in particular, the filters would be washed one time with 1 x SSPE. At this point procedures discussed later in the application apply.

The probe is designed and selected to provide a complementary sequence for binding to the target DNA or RNA sequences. Short-chain probes are preferred, such as those containing from 10 to 50

nucleotide units, and a single complementary sequence to a characterizing sequence of the target DNA or RNA. Longer probes such as probes containing in excess of 200 nucleic acid units may contain a plurality of complementary sequences. Such long probes may be used to obtain some of the advantages of the present invention. Probes prepared by nucleotide synthesis which contain from 18 to 30 nucleic acid units are particularly desirable. Methods for synthesizing and labelling such short probes are described in the published PCT Application WO 84/03285.

As in prior practice, any one of a number of different kinds of reporter groups or labels may be attached to the probes. Radioactive labels which may be employed include ^{32}P , ^3H , ^{14}C labelled probes. However, labels generating some fluorescent or colorimetric response are preferred, as described in the published European Patent Application 0 070 687. In particular, it is preferred to employ fluorescent labels. Fluorescent labelling is described in PCT Publication No. WO 84/03285. The fluorescent labels may include fluorescein, Texas Red, Lucifer Yellow, pyrene, lanthanide complexes, etc. The preferred probe is one containing a single fluorescent label attached to a complementary sequence of 18 to 30 nucleic acid units. However, the use of two or more fluorophore labels per probe is not excluded, particularly in cases where fluorophore combinations effect a condition for improved detection of fluorescent signal.

The novel steps of the present invention begin following the hybridization and removal of the unbound probe from the initial support matrix (nitrocellulose, etc.). In the event that the concentration of the hybridized probe is sufficient to provide appropriate selectivity and sensitive detection, the method of the present invention need not be employed.

The first step in utilization of the new methodology involves the dehybridization of the labelled probe from the initial support matrix. Dehybridization may be effected by known denaturation procedures, which cause double-stranded nucleic acids to separate into single-stranded form. For example, dehybridization may be carried out by heating the support matrix-target DNA/hybridized labelled probe in a small amount of water or dilute buffer (0.1 to 1X SSPE or SSC) to a temperature sufficient to produce denaturation. In general, effective denaturing temperatures (also called melting temperatures) are in the range from 60° to 90°C. The required temperature depends on the length and "GC" (Guanine & Cytosine) content of the probe. In an alternative procedure, an appropriate buffered solution containing a denaturing agent such as urea or formamide can be used to cause dehybridization of the probe. For example, a 6 molar urea (0.1X - 1X SSPE or SSC) solution is effective in dehybridizing the probe at room temperature (20°-30°C). A variety of general procedures and agents exist for causing dehybridization of the probe from the support. In general, selection is based on the ease and speed of the procedure and the important requirement that it not interfere or inhibit the second stage of the process, i.e., where the labelled probe is adsorbed and concentrated onto a suitable material for detection.

In an alternative procedure which is not known to have been previously employed, the dehybridization of the labelled probe is effected by means of a "displacer probe." The displacer comprises a single-stranded oligonucleotide sequence including the same sequence as the labelled probe with an additional sequence on one side, the additional sequence being complementary to a corresponding sequence on the target strand to which the probe is bound. The additional sequence should be of sufficient length to enable the displacer to preferentially bind to the target sequence near the probe and then efficiently displace the probe. For example, the additional sequences may contain from 10 to 20 nucleic acid units and may be provided on either side of the labelled probe sequence. By way of more specific example, where the hybridized labelled-probe contains from 18 to 30 nucleic acid units, the displacer may correspondingly contain from about 28 to 50 nucleic acid units. Longer displacer probes can also be employed.

In using the displacer, it may be applied in aqueous buffered solution to the hybridized probe. For example, the supported hybridized probe may be immersed in a solution of the displacer. The typical concentrations will range from about 1×10^{-13} to 1×10^{-11} moles of the displacer per 50 microliters of solution (0.1 to 5.0 x SSPE or SSC). The displacement may be carried out at room temperatures (viz. 20-30°C) in less than 30 minutes. In general, the conditions employed are optimized to allow efficient hybridization of the displacer probe/dehybridization of the labelled probe, and as minimal leaching of non-specific bound labelled probe as possible.

This displacement procedure has a further important advantage in that it increases selectivity. During hybridization some of the labelled probe may have become bound to components of the sample other than the target sequence. Where such extraneous binding has occurred, a false positive indication can be obtained if the non-specifically bound labelled probe is detected in the sample. By using the displacer, however, displacement occurs for only that portion of the probe specifically hybridized to the target sequence. Thus, the displacer probe gives double specificity to the assay procedure.

With the preferred procedures described above, the bound labelled probe can be dehybridized from the support material and obtained in the form of an aqueous solution. Direct detection of the labelled probe in

the aqueous solution is usually not feasible because the concentration of the labelled probe is extremely low. However, by concentrating and focusing the dehybridized probe, it can be effectively detected at concentrations down to 10^{-8} moles and below.

In one preferred procedure, the dehybridized probe is adsorbed from the solution onto an anion exchange material which provides positively charged groups for binding the probe. Weak basic anion exchange materials containing substituted amino groups may be employed as well as strongly basic anion exchange materials with quaternary amino groups.

It is presently believed, however, that quaternary amino type materials are the most desirable.

Anion exchange materials suitable for use in the method of this invention are available commercially. Such materials can be prepared from substrates comprising cross-linked dextran or glass. Such materials are usually in the form of porous beads or granules of ultra-fine particle sizes (viz. 50-200 micron diameter). Anion exchange materials based on cross-linked dextran are available from Pharmacia (Uppsala, Sweden). These products are in the form of small, generally spherical beads with varying degrees of porosity depending on the extent of cross-linking of the dextran. Quaternary amino beads are sold by Pharmacia under the trademark "QAE-Sephadex." The quaternary groups are attached to the dextran substrate as quaternary amino-ethyl groups. Other amine-type beads are also supplied by Pharmacia under the trademark "DEAE-Sephadex." These products contain diethylaminoethyl groups attached to the dextran.

While either DEAE-Sephadex or QAE-Sephadex can be used with the method of the present invention, it is presently believed that "QAE-Sephadex" is preferred. Further, a presently preferred adsorbent is "QAE-Sephadex A-25," which is more highly cross-linked dextran substrate than "QAE-Sephadex A-50." Greater cross-linking restricts swelling of the beads and this is believed advantageous since an essentially surface adsorption is desired for purposes of the present invention. However, some penetration of the probe into the beads is not objectionable. Typically, the Sephadex beads may have average sizes of around 100 microns.

Where an adsorbent is employed for the dehybridized probe as described above, it is desirable that the adsorbent be formed of a substance which does not interfere with the fluorescent analysis. Fortunately, it has been found that both anion exchange dextran beads and glass beads have excellent background properties with reference to fluorescent analysis. Neither the dextran nor the glass interfered with the analysis. These adsorbents have very low background fluorescence and produce minimal light scattering effects.

Porous glass beads are also available which have anion exchange properties, being either weak ion exchangers with attached amine groups or strong ion exchangers with attached quaternary amino groups. For example, such products can be obtained from Pierce Chemical Company, Rockford, Illinois. The size of the beads may range from about 50 to 200 microns. This type of product is sometimes referred to as "controlled pore glass." Diethylaminoethyl (DEAE) controlled pore glass is available from Pierce Chemical Company as well as quaternary aminoethyl (QAE) controlled pore glass. One specific product which may be employed is sold by Pierce as Product No. 23514 CPG/QAE Glycophase Pore Glass.

For use in conjunction with fluorescent detection equipment, such as an epifluorescent microscope, it may be desirable to prepare glass slides, glass tips, or fiber optic tips which provide a surface that is suitable for binding the dehybridized probe. Amine groups can be attached to the silica glass by a silanization process. See, for example, Lowe and Dean, "Affinity Chromatography, pp. 217 and 256 (1974), John Wiley & Sons, New York. As initially attached, the nitrogen groups will be in the form of primary amine groups. By well-known procedures, such primary amine groups can be easily converted to secondary, tertiary, or quaternary amino groups.

To achieve effective concentration of the dehybridized probe, it is preferred to use minimal amounts of the adsorbent. For example, employing Sephadex QAE-A25 beads (~100 micron diameter), it will usually be sufficient to use from 10 to 50 of the beads to adsorb labelled probe from 20-50 μ l of solution (probe originally hybridized to a two millimeter diameter sample dot on a nitrocellulose filter). Usually about 50% or more of the labelled probe can be adsorbed onto the beads in less than one hour. In the case of using 10 beads, about a 40 fold concentration of labeled probe is achieved in moving it from the 2 millimeter nitrocellulose sample dot (area $\sim 1 \times 10^5$ square microns). With the probe in a 50 μ l volume (5×10^9 cubic microns) solution adsorption to 10 beads (1×10^7 cubic microns) produces a 5000 fold concentration effect in terms of volume. In general, the smaller the number of beads or surface area the probe is finally adsorbed to the better. However, for any given number of beads there is a limitation of the probe transfer rate to the beads as the volume of the solution increases. Since the adsorbent materials are ion-exchangers, efficiency for binding probes will decrease as ionic strength of the solutions increases. Therefore solutions should be kept as diluted as possible (≤ 0.4 M NaCl and ≤ 0.02 M for phosphate and citrate buffers). In general, extremes of temperature or pH should be avoided with regard to adsorption of

labelled probe to beads. Temperatures higher than $\sim 40^{\circ}\text{C}$ and pH's lower than 5 should be avoided in particular.


It is within the scope of this invention to employ other methods for concentrating and/or increasing the transfer rate of the dehybridized probe to adsorbent materials. For example, the solution containing the probe may be concentrated by evaporation, flocculation or precipitation. These and similar procedures can result in concentration and focusing of the dehybridized probe permitting improved fluorescent analysis. With regard to increasing or improving the transfer rate of labelled probe in solution to the adsorbent material, general procedures include agitation and mixing of the adsorbent beads with the labelled probe solution or the flowing of the labelled probe solution through an aggregation of the adsorbent material in a small orifice. Also, it is possible to electrofocus the labelled probe onto a small area using microelectrophoresis.

After the dehybridized probe has been adsorbed and concentrated, such as on the adsorbent beads described above, it can be examined by standard procedures for detecting the particular label or reporter group associated with the probe. Where the probe, as preferred, is labelled with a fluorescent group, the beads may be examined with an epifluorescent microscope or other fluorescent analysis system designed for ultrasensitive detection. Usually the analysis and quantitation takes only a matter of minutes for a given sample.

The method of this invention can be further understood in relation to the following examples.

EXAMPLE I

The following example demonstrates the advantage for detecting fluorescent probes concentrated onto suitable materials, like Sephadex QAE A-25 beads, over observing the same probe bound to nitrocellulose filter material. The fluorescent probe used in these experiments was a 22-mer sequence specific for a genomic target sequence present in Herpes Simplex Virus (HSV). The probe was labeled with a Texas Red fluorophore group having an excitation and emission maxima at 596nm and 620nm respectively. The labeled probe was designated TR-HSV, and the sequence and position of the label are shown below:

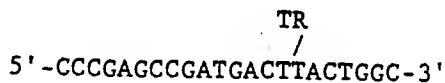

The first experiment involved the epifluorescent analysis of the TR-HSV probe bound to nitrocellulose filters. A 10-fold dilution series of TR-HSV probe ranging from 2×10^{-12} to 2×10^{-18} mole per microliter of 1X SSPE buffer was prepared. One microliter of each dilution was adsorbed onto nitrocellulose filter material forming a 2-3 millimeter diameter spot. A control spot of $1\mu\text{L}$ of 1 x SSPE (no probe) was also prepared. The spots were now analyzed by using a photon counting epifluorescent microscope filtered for Texas Red fluorescence (Ex $\sim 550\text{nm}$; Em $\sim 630\text{nm}$). The fluorescent signal from each sample spot was counted (EG&G Photon Counting System) for one second (each sample spot was counted twice). Results for the analysis are shown in Table I. Since only 1/10 the area of each spot was analyzed (10x microscope field), samples 1 through 7 range from 2×10^{-13} to 2×10^{-19} mole in actual TR-HSV probe concentration.

TABLE I

TR-HSV Probe/Nitrocellulose Filter

Sample	TR-HSV (mole)*	Avg. Counts (1 second)
1	2×10^{-13}	1.5×10^6
2	2×10^{-14}	157,429
3	2×10^{-15}	27,918
4	2×10^{-16}	26,246
5	2×10^{-17}	21,155
6	2×10^{-18}	22,100
7	2×10^{-19}	23,397
Control	---	22,399

*Concentration in area observed

The results in Table I show that fluorescence from the TR-HSV probe can be detected down to about 2×10^{-16} mole (sample 4) on nitrocellulose filter material.

The second experiment involved the epifluorescent analysis of the TR-HSV probe bound to QAE-Sephadex A-25 beads. A 10-fold dilution series of TR-HSV probe ranging from 4×10^{-14} to 4×10^{-20} mole per 10 microliters of 1 X SSPE buffer was prepared. Approximately 1000 QAE-Sephadex A-25 beads (previously swollen in 1 X SSPE buffer) were added to each of the sample dilutions. A control sample of beads with no TR-HSV probe was also prepared. The samples were allowed to stand for about one hour with intermittent gentle agitation. Samples were then analyzed using the photon counting epifluorescent microscope. In this experiment, the fluorescence from individual beads in each sample were counted. The results in Table 2 are the average counts (one second) from two beads in each sample counted twice. The results are given in terms of TR-HSV concentration per bead, i.e., it is assumed most of the probe has adsorbed evenly to all the beads in the sample.

TABLE 2

TR-HSV Probe/QAE-Sephadex A-25

Sample	TR-HSV (mole)/bead	Avg. Counts (1 second)
1	4×10^{-14}	9×10^6
2	4×10^{-15}	645,867
3	4×10^{-16}	43,229
4	4×10^{-17}	15,614
5	4×10^{-18}	13,581
6	4×10^{-19}	5,429
7	4×10^{-20}	4,782
Control	---	4,320

The results in Table 2 show that fluorescence from the TR-HSV probe at 4×10^{-20} mole per bead still has counts significantly higher than the control beads (no probe).

The third experiment is basically a complete repeat of the second experiment. The results are given in Table 3.

TABLE 3

TR-HSV Probe/QAE Sephadex A-25

Sample	TR-HSV (mole)/bead	Avg. Counts (1 second)
1	4×10^{-14}	1×10^7
2	4×10^{-15}	1×10^6
3	4×10^{-16}	85,939
4	4×10^{-17}	17,022
5	4×10^{-18}	13,627
6	4×10^{-19}	10,182
7	4×10^{-20}	8,082
Control	---	7,341

The results in Table 3 show again the fluorescence from the TR-HSV probe at 4×10^{-20} mole per bead has counts significantly (~10%) higher than the control beads. The higher counts for control beads in experiment 3 (7,341) compared to experiment 2 (4,320) are due to the higher excitation intensity used in this experiment.

In summary, TR-HSV probe can be detected down to about 4×10^{-20} mole on a single bead of QAE-Sephadex A-25 by epifluorescence analysis. On nitrocellulose filter material it can be detected down to only about 2×10^{-16} mole.

EXAMPLE II

In further experiments similar to those described in Example 1 it was demonstrated that another type of fluorescent probe (Fluorescein labelled) can be concentrated and more favorably detected on QAE-Sephadex A-25 beads. Also, it is shown that other supports, such as QAE Controlled Pore Glass Beads, work favorably.

The first experiment involves the epifluorescent analysis of a fluorescein labelled HSV probe (F-HSV) bound to Zeta Bind (nylon) filter material. The F-HSV probe is exactly the same as the TR-HSV probe except for the different fluorophore label. The fluorescein label has its excitation and emission maximum at ~490nm and ~520nm respectively. The experiment with F-HSV probe adsorbed to Zeta Bind was carried out exactly like the first experiment in Example I. However, the Photon counting epifluorescence microscope was filtered for fluorescein fluorescence (Ex ~ 450 nm, Em ~ 520 nm). Results for the analysis are shown in Table 4.

TABLE 4

F-HSV Probe/Zeta Bind Filter

Sample	F-HSV (mole)*	Avg. Counts (1 second)
1	6×10^{-13}	843,568
2	6×10^{-14}	200,248
3	6×10^{-15}	55,914
4	6×10^{-16}	45,587
Control	---	47,675

*Concentration in area observed.

The results from Table 4 show that fluorescence at 6×10^{-15} mole of probe is just distinguishable from the control sample.

The second experiment involves the epifluorescent analysis of the F-HSV probe bound to QAE-Sephadex A-25 beads. This experiment was carried out in a manner similar to the second experiment in Example I. Results for the analysis are shown in Table 5.

TABLE 5

F-HSV Probe/QAE-Sephadex Beads

Sample	F-HSV (mole)/bead	Avg. Counts (1 second)
1	7×10^{-14}	10^7
2	7×10^{-15}	2×10^6
3	7×10^{-16}	136,736
4	7×10^{-17}	56,924
5	7×10^{-18}	31,401
6	7×10^{-19}	32,194
Control	---	27,507

The results from Table 5 shown that fluorescence from the F-HSV probe is detectable at the 7×10^{-18} mole per bead level.

As was the case with the TR-HSV probe (Example I), the F-HSV probe is significantly more detectable on QAE-Sephadex A-25 beads than on Zeta Bind filter material.

The third experiment involves the epifluorescent analysis of TR-HSV probe bound to QAE-Controlled Pore Glass beads. This experiment was carried out similar to the others described above. Results for the analysis are shown in Table 6.

TABLE 6

TR-HSV Probe/QAE-Controlled Pore Glass Beads

Sample	TR-HSV Probe/Bead	Avg. Counts (1 second)
1	4×10^{-14}	10^7
2	4×10^{-15}	2×10^6
3	4×10^{-16}	179,131
4	4×10^{-17}	51,797
5	4×10^{-18}	20,454
6	4×10^{-19}	16,897
7	4×10^{-20}	15,497
Control	---	11,930

The results from Table 6 show that TR-HSV probe is detectable to the 4×10^{20} mole level on individual QAE-Controlled Pore Glass beads.

EXAMPLE III

The experiment described below demonstrates that fluorescent probe actually hybridized to a target DNA sequence immobilized on a polystyrene bead support material (not suitable for fluorescent analysis) can be dehybridized and concentrated onto QAE-Sephadex A-25 beads for epifluorescent analysis. In this experiment the TR-HSV probe used was also labelled with ^{32}P Phosphorous radioisotope (^{32}P). Thus a comparison could be made between the fluorescent and radioisotope level of detection.

Using one milligram samples of polystyrene beads (~1-2 micron diameter, substituted with a complementary target sequence), a 10-fold dilution series of ^{32}P -TR-HSV probe ranging from 1×10^{-13} to 1×10^{-18} mole in 50 μl of 5X SSPE hybridization buffer were prepared (control sample contained no probe). The samples were hybridized for 30 minutes at 40°C , then centrifuged and supernatant removed. The samples were then washed three times with 5X SSPE. After the final wash and removal of supernatant, 10 μl of 6M Urea (1X SSPE) was added to each of the samples in order to dehybridize the probe from the support. Dehybridization was carried out for about 30 minutes at room temperature. At this point the polystyrene beads were separated from the solution containing dehybridized ^{32}P -TR-HSV probe by centrifugation. To this solution between 20 to 30 QAE-Sephadex A-25 beads were added. The dehybridized probe was allowed to adsorb to the QAE-Sephadex beads for about 30 minutes. At this point the QAE-Sephadex A-25 beads were separated from the solution. The beads (all 20 to 30) were first counted on a scintillation counter (one minute counts) to determine level of probe finally transferred to the QAE-Sephadex A-25 beads by radioisotope method. After scintillation counting individual beads were counted on the photon counting epifluorescent microscope to determine sensitivity by fluorescence. The results for each type of analysis are shown in Table 7.

TABLE 7

³²P-Tr-HSV Transfer to QAE-Sephadex

Sample	Total Conc. ³² P-Tr-HSV Probe	³² P Counts	³² P-Tr-HSV Conc. per bead	Fluorescent Counts (1 sec.)
1	$\sim 6 \times 10^{-15}$	7460	$\sim 3 \times 10^{-16}$	104,572
2	$\sim 6 \times 10^{-16}$	951	$\sim 3 \times 10^{-17}$	28,847
3	$\sim 6 \times 10^{-17}$	159	$\sim 3 \times 10^{-18}$	13,704
4	$\sim 6 \times 10^{-18*}$	43	$\sim 3 \times 10^{-19}$	13,170
5	$\sim 6 \times 10^{-19*}$	40	$\sim 3 \times 10^{-20}$	10,227
Control	---	44	---	8,900

*extrapolated

The results in Table 7 show that about 6×10^{-17} mole of ³²P-TR-HSV probe transferred to QAE-Sephadex beads could be detected by radioisotope counting (³²P, 1 minute count) before control or background level is reached (~ 44 counts). Fluorescent counts were found to be above the control through the whole sample series 1 \rightarrow 5. When considered in terms of individual beads, on the order of 10^{-20} mole of ³²P-TR-HSV probe was detectable.

EXAMPLE IV

The following experiment demonstrates the usefulness of a so-called "displacer probe" in effecting the dehybridization of a fluorescence labelled probe hybridized to target DNA immobilized on a support material.

In this experiment ³²P-TR-HSV probe was hybridized to plasmid target DNA (pHSV106) immobilized onto nitrocellulose filter material. The pHSV106 plasmid DNA was immobilized onto the nitrocellulose filter using the Hybrid-Dot apparatus and experimental procedure described earlier in this patent application. The concentrations of pHSV106 plasmid DNA in each "dot" (~ 2 millimeter diameter spot) were 1) 100ng; 2) 10 ng, 3) 1 ng, 4) 100 pg, and 5) control (non-target DNA). ³²P-TR-HSV probe was hybridized to the pHSV106 target dots by the procedures discussed earlier in this application. After hybridization and appropriate washing, the level of ³²P-TR-HSV probe hybridized to each "dot" was measured by a scintillation counter. The results are given in Table 8. After counting, the dehybridization procedure was carried out using a specific displacer probe for the ³²P-TR-HSV probe. The sequence of the displacer probe is given below:

5'-CCCCAGCACCTGCCAGTAAGTCATCGGCTCGGG-3'

The portion of the displacer probe underlined is identical to that of the ³²P-TR-HSV probe. Dehybridization of ³²P-TR-HSV/pHSV106 target DNA nitrocellulose "dots" was carried out by treating each one of the sample "dots" with 50 μ l of 1X SSC buffer solution containing approximately 100ng of the displacer probe for 30 minutes at $\sim 25^{\circ}\text{C}$. The nitrocellulose filter "dots" were then separated from the solution and each solution and "dot" sample was counted on the scintillation counter. The results are given in Table 8.

TABLE 8
Displacer Probe Experiment

5

	Sample (pHSV106)	1st Count (1 min.) after Hybridization & Wash	2nd counts (1 min) After Displacer Probe	
			Filter "Dots"	Solution
10	1 100ng	39,672	5,969	32,353
	2 10ng	5,805	1,347	4,296
	3 1ng	1,702	1,113	754
15	4 100pg	1,543	1,109	706
	5 None	1,296	1,210	198

20 The results in Table 8 show that the "displacer probe" treatment has successfully dehybridized the ³²P-TR-
HSV probe from the pHSV106 target/nitrocellulose filter "dots."

Claims

25 1. The hybridization assay method for detecting a specific polynucleotide target sequence wherein a nucleic acid containing target sample in single stranded form is affixed to a support, contacted thereon under hybridizing conditions with a reporter group-labeled single-stranded polynucleotide probe having a sequence complementary to the target sequence, at least part of said probe being bound to said affixed sample when the target is present, said probe hybridizing with said target sequence, then removing the unbound portion of said probe from said support, and detecting the presence of said probe by means of its reporter group, wherein the sensitivity of the assay is increased by employing the additional steps

30 comprising:
 (a) after removal of the unbound portion of the probe, dehybridizing the bound portion;
 (b) separating the dehybridized portion from the support and the residual sample; and
 35 (c) concentrating the dehybridized separated probe before detecting its presence by means of said reporter group.

2. The method of claim 1 in which the dehybridized portion of the probe is concentrated by adsorption from an aqueous solution onto a basic anion exchange material binding said probe.

3. The method of claim 2 in which the anion exchanger functional groups are quaternary amine groups.

40 4. The method of claim 3 in which said anion exchange functional groups are provided on a substrate selected from cross-linked dextran or glass.

5. The method of claim 1 in which the bound labelled probe portion is dehybridized by contacting the supported sample with a displacer probe comprising a single stranded polynucleotide sequence including the same sequence as the labelled probe with an additional sequence on at least one side thereof complementary to the corresponding sequence on the target sequence to which the labelled probe is bound, said additional sequence being of sufficient length to cause said displacer to preferentially bind to said target sequence and subsequently displace said labelled probe therefrom.

6. The method of claim 1 in which said probe comprises a sequence of 10 to 50 nucleotide units and has a reporter group attached thereto, and said probe is concentrated on said anion exchange material before detecting its presence by means of said reporter groups.

7. The method of claim 6 in which said anion exchange material comprises an adsorbent substrate selected from cross-linked dextran or glass with substituted amino groups attached thereto, and is in the form of generally spherical beads.

8. The method of claim 6 in which said probe is labelled with a fluorescent reporter group.

55 9. The method of claim 6 in which the bound probe portion is dehybridized by contacting the supported sample with a displacer comprising a single-stranded polynucleotide sequence including the same sequence as the probe with an additional sequence on at least one side thereof complementary to the

corresponding sequence on the same side of the target sequence to which the probe is bound, said additional sequence being of sufficient length to cause said displacer to preferentially bind to said target sequence and displace said probe therefrom.

10. The method of claim 6 in which the bound probe portion is dehybridized by contacting the supported sample with a displacer comprising a single stranded polynucleotide sequence including the same sequence as the probe with an additional sequence on at least one side thereof complementary to the corresponding sequence on the same side of the target sequence to which the probe is bound, said additional sequence being of sufficient length to cause said displacer to preferentially bind to said target sequence and displace said probe therefrom.

10

15

20

25

30

35

40

45

50

55



European Patent
Office

EUROPEAN SEARCH REPORT

Application number

EP 86 11 8191

DOCUMENTS CONSIDERED TO BE RELEVANT			
Category	Citation of document with indication, where appropriate, of relevant passages	Relevant to claim	CLASSIFICATION OF THE APPLICATION (Int. Cl. 4)
Y	EP-A-0 142 299 (FUJIREBIO K.K.) * abstract *	1-8	C 12 Q 1/68
Y,D	US-A-4 358 535 (S. FALKOW et al.) * column 1, line 64 - column 2, line 61; column 3, line 25 - column 4, line 4 *	1-8	
A	EP-A-0 167 238 (ALLIED CORP.) * claims 1-3 *	1,5	
A	EP-A-0 171 072 (AKADEMIE DER WISSENSCHAFTEN DER DDR) * pages 3-4, complete; claims 1-4 *	1-3	
A	ANALYTICAL CHEMISTRY, vol. 58, no. 13, November 1986, pages 1306-1310, Washington DC, US; S.A.B. "DNA Probes"		TECHNICAL FIELDS SEARCHED (Int. Cl. 4) C 12 Q 1/00 G 01 N 33/00
The present search report has been drawn up for all claims			
Place of search BERLIN		Date of completion of the search 22-07-1987	Examiner GREEN C.H.
CATEGORY OF CITED DOCUMENTS X : particularly relevant if taken alone Y : particularly relevant if combined with another document of the same category A : technological background O : non-written disclosure P : intermediate document T : theory or principle underlying the invention E : earlier patent document, but published on, or after the filing date D : document cited in the application L : document cited for other reasons & : member of the same patent family, corresponding document			

(19)



Europäisches Patentamt
European Patent Office
Office européen des brevets



(11) Publication number:

0 269 764 B1

(12)

EUROPEAN PATENT SPECIFICATION

(45) Date of publication of patent specification: 28.08.91 (51) Int. Cl.⁵: C12Q 1/68

(21) Application number: 86118191.5

(22) Date of filing: 01.12.86

(54) Method for increasing the sensitivity of nucleic acid hybridization assays.

(43) Date of publication of application:
08.06.88 Bulletin 88/23

(45) Publication of the grant of the patent:
28.08.91 Bulletin 91/35

(84) Designated Contracting States:
AT BE CH DE ES FR GB GR IT LI LU NL SE

(56) References cited:
EP-A- 0 142 299
EP-A- 0 167 238
EP-A- 0 171 072
US-A- 4 358 535

ANALYTICAL CHEMISTRY, vol. 58, no. 13, November 1986, pages 1306-1310, Washington DC, US; S.A.B. "DNA Probes"

(73) Proprietor: MOLECULAR BIOSYSTEMS, INC.
11180 Roselle Street, Suite A
San Diego, CA 92121(US)

(72) Inventor: Heller, Michael J.
1357 Acton Avenue
Poway, CA 92064(US)

(74) Representative: Patentanwälte Grünecker,
Kinkeldey, Stockmair & Partner
Maximilianstrasse 58
W-8000 München 22(DE)

EP 0 269 764 B1

Note: Within nine months from the publication of the mention of the grant of the European patent, any person may give notice to the European Patent Office of opposition to the European patent granted. Notice of opposition shall be filed in a written reasoned statement. It shall not be deemed to have been filed until the opposition fee has been paid (Art. 99(1) European patent convention).

Description

The field of this invention involves a hybridization assay system for detecting target polynucleotide sequences in deoxyribonucleic acid (DNA) and/or ribonucleic acid (RNA) samples. More specifically, the invention is concerned with a method for increasing the sensitivity and where necessary the selectivity of the hybridization assay using labelled oligonucleotide probes to detect specific target sequences.

Hybridization assays may be employed for the detection and identification of DNA or RNA sequences. Published methods as used particularly in recombinant DNA research are described in *Methods in Enzymology*, Vol. 68, pp. 379-469 (1979); and Vol. 65, Part 1, pp. 468-478 (1968). One such method involving a preliminary separation of the nucleic acid fragments by electrophoresis is known as the "Southern Blot Filter Hybridization Method." See Southern, *J. Mol. Biol.* (1975) 98:503. Hybridization probe methods for detecting pathogens are described in U.S. Patent 4,358,539.

There has been a recognized need for improving hybridization assays. Instead of radiolabelled probes, which can require several days for development of the X-ray film, detection by means of luminescent-labelled probes are claimed to have advantages. (See published European Patent Application 0 070 687.) A procedure for obtaining increased sensitivity from luminescent-labelled probes is disclosed in published European Patent Application 0 070 685. Two probes are employed, having respectively a chemiluminescent catalyst and an absorber/emitter (fluorophore) moiety. Upon hybridization to the appropriate polynucleotide target sequence, the light signal from the probes is shifted to a higher wavelength. This wavelength shifting leads to an improved and more sensitive assay system.

As an alternative to employing long double-stranded probes resulting from recombinant DNA cloning procedures, short synthetically produced single-stranded probes have been described. (See, for example, PCT Publication Number WO 84/03285). Such short probes, which normally contain from 18 to 40 nucleotides, are available from Molecular Biosystems, Inc., San Diego, California. Typically, such short probes can be labelled with biotin, enzymes, and fluorescent reporter groups. These probes can be constructed complementary to a unique sequence of the target DNA or RNA fragments. Furthermore, they are in stable single-stranded form and renaturation cannot interfere with the assay. Short labelled DNA probes as described above, in principle, should be useful in developing rapid and convenient non-radioisotope DNA hybridization assays. With existing procedures, concentrations of the hybridized labelled probe can be detected down to about 10^{-16} to 10^{-17} moles under ideal conditions. For most of the potentially important clinical applications of DNA probes (infectious diseases, latent virus, genetic disorders, etc.), it will be necessary to have sensitivity levels of at least 10^{-18} moles and lower. With present techniques the inherent sensitivities possible with non-radioisotope labels such as fluorophores and enzymes are still limited by a wide variety of problems. In the case of fluorescent labels, background fluorescence and light scattering from the hybridization support materials (nitrocellulose and nylon filters, polystyrene beads, etc.) lead to significant loss of inherent sensitivity. In the case of enzyme labels, non-specific binding of the labeled probe and interfering stains and colors from sample material can limit ultimate sensitivity of this label. Other problems superimposed upon those mentioned above lead to further losses of sensitivity and limit the practicality of these systems for clinical diagnostic purposes.

SUMMARY OF INVENTION

The present invention provides a method for ultrasensitive detection of reporter group-labelled DNA probes used in hybridization assay systems. The method is particularly applicable to fluorescently-labelled short probes, but should also provide improvements for enzyme-labelled probes. In the case of fluorescent probes, greater sensitivities can be obtained than with any prior hybridization assay system, viz. as little as 10^{-18} moles of probe and below can be detected.

The method of this invention utilizes a probe concentration or refocusing procedure. After the basic hybridization procedure has been carried out and the support matrix (nitrocellulose, for example) has been washed free of unbound probe, instead of attempting to detect the presence of bound probe directly as in conventional systems, the probe is subjected to a dehybridization procedure. Dehybridization can be carried out by treatment with denaturing agents (heat, urea, etc.) which will release the bound probe from the support matrix and liberate it into the aqueous solution. The bound probe can also be liberated in a much more selective fashion using a second DNA probe which is called a "displacer probe." In either case the liberated probe, in solution, is now concentrated or refocused by addition of a small amount of an appropriate material with high affinity for the probe.

In one preferred procedure, after dehybridization from the support matrix (nitrocellulose or nylon filters), the fluorescently labelled probe is adsorbed onto a small amount of adsorbent beads or particles having

high affinity for oligonucleotides. For example, the adsorbent may comprise beads of porous, crosslinked dextran or glass particles substituted with positively charged groups. Such adsorbents are available commercially and are improved substrates for fluorescent probe detection in particular. For example, anion exchange materials such as diethylaminoethyl (DEAE) or quaternary aminoethyl (QAE) dextran and glass beads have extremely high affinity for labelled oligonucleotide probes. This high affinity is due to the electrostatic or ion exchange interaction between the polyanionic nucleotide and the high positive charge density on the beads. These beads have been found to provide extremely low background fluorescence and significantly less light scattering effects than nitrocellulose or nylon filter material. Most importantly, the dehybridization and adsorption onto a small amount of adsorbent material concentrates (10-5000 fold) the fluorescent probe onto a substrate much more favorable for fluorescent (epifluorescent) detection. In summary, the fluorescent signal is now highly concentrated, background is significantly reduced, and interferences from the hybridization support material (nitrocellulose filters, other sample components on the filters, etc.) have been eliminated.

While the specific procedures discussed above are preferred, within the broader context of the improved method of this invention other concentration or focusing procedures can be employed. Moreover, the probes can be used which are labelled with other types of reporter groups (i.e., enzymes or luminescent groups). With such modifications, improvements in detection sensitivity can also be obtained in accordance with the principles of this invention.

DETAILED DESCRIPTION

In practicing the method of this invention, the DNA or RNA test samples containing target sequences may be prepared by any one of a number of known procedures and attached to suitable immobilization support matrix. For example, such procedures are described in Methods in Enzymology, Vol. 68, pp. 379-469 (1979); and Vol. 65, Part 1, pp. 468-478 (1980), U.S. Patent 4,358,539; and in published European Patent Applications Nos. 0 070 685 and 0 070 687. Any of the usual supports can be employed, such as nitrocellulose filters, Zetabind® (nylon) filters, polystyrene beads, etc.

The following is a general description of the basic procedures for carrying out the hybridization of labelled short DNA probes to samples of target DNA immobilized on nitrocellulose and Zetabind® filters. With minor variations the procedure is applicable for hybridizing radioisotope (³²P), biotin, enzyme, lumiphore, and fluorophore labelled short probes.

Buffers and Solutions

All buffers should be made with double-distilled water. Concentrated stock solutions should be diluted with double-distilled water. Filter, where indicated, through 0.45µ nitrocellulose.

20 SSPE BUFFER. Dissolve 210 g NaCl, 27.6 g NaH₂PO₄ · H₂O, 7.4 g Na₂ EDTA, and approximately 8.3 ml of saturated NaOH (to pH 7.0). Adjust volume to one liter, filter, autoclave and store at room temperature.

20 SSC BUFFER. Dissolve 175.3 g NaCl and 88.2 g sodium citrate in 800 ml water. Adjust pH to 7.0 with a few drops of concentrated NaOH. Adjust volume to one liter, filter, autoclave, and store at room temperature.

0.3 M HAcOH. Store at room temperature.

1 M and 2 M Ammonium Acetate. Filter and store at room temperature.

4 x BP. Dissolve 2 g each bovine serum albumin, crystalline Pentex Fraction V., and polyvinylpyrrolidone, average mol. wt. = 40,000 in 100 ml water. Filter and store at -20 °C in small aliquots.

10% (W/V) Sodium Dodecyl Sulfate (SDS). Store at room temperature.

Other Materials:

Nitrocellulose filters (Schleicher & Schuell, Keene, NH, #BA85) or Zetabind membranes (AMF, Meriden, CT, #NM511-01-04-5SP).

Whatman 3MM chromatography paper.

DEAE-cellulose. Whatman DE52, precycle according to the manufacturer.

Dot blot filtration apparatus (HYBRIDOT™, BRL, Bethesda, MD; Minifold II, S & S, Keene, NH; BIODOT™, Bio-Rad Laboratories, Richmond, CA.)

Water baths.

Vacuum pump.

high affinity for oligonucleotides. For example, the adsorbent may comprise beads of porous, crosslinked dextran or glass particles substituted with positively charged groups. Such adsorbents are available commercially and are improved substrates for fluorescent probe detection in particular. For example, anion exchange materials such as diethylaminoethyl (DEAE) or quaternary aminoethyl (QAE) dextran and glass beads have extremely high affinity for labelled oligonucleotide probes. This high affinity is due to the electrostatic or ion exchange interaction between the polyanionic nucleotide and the high positive charge density on the beads. These beads have been found to provide extremely low background fluorescence and significantly less light scattering effects than nitrocellulose or nylon filter material. Most importantly, the dehybridization and adsorption onto a small amount of adsorbent material concentrates (10-5000 fold) the fluorescent probe onto a substrate much more favorable for fluorescent (epifluorescent) detection. In summary, the fluorescent signal is now highly concentrated, background is significantly reduced, and interferences from the hybridization support material (nitrocellulose filters, other sample components on the filters, etc.) have been eliminated.

While the specific procedures discussed above are preferred, within the broader context of the improved method of this invention other concentration or focusing procedures can be employed. Moreover, the probes can be used which are labelled with other types of reporter groups (i.e., enzymes or luminescent groups). With such modifications, improvements in detection sensitivity can also be obtained in accordance with the principles of this invention.

20 DETAILED DESCRIPTION

In practicing the method of this invention, the DNA or RNA test samples containing target sequences may be prepared by any one of a number of known procedures and attached to suitable immobilization support matrix. For example, such procedures are described in *Methods in Enzymology*, Vol. 68, pp. 379-469 (1979); and Vol. 65, Part 1, pp. 468-478 (1980), U.S. Patent 4,358,539; and in published European Patent Applications Nos. 0 070 685 and 0 070 687. Any of the usual supports can be employed, such as nitrocellulose filters, Zetabind® (nylon) filters, polystyrene beads, etc.

The following is a general description of the basic procedures for carrying out the hybridization of labelled short DNA probes to samples of target DNA immobilized on nitrocellulose and Zetabind® filters. With minor variations the procedure is applicable for hybridizing radioisotope (32p), biotin, enzyme, lumiphore, and fluorophore labelled short probes.

Buffers and Solutions

All buffers should be made with double-distilled water. Concentrated stock solutions should be diluted with double-distilled water. Filter, where indicated, through 0.45µ nitrocellulose.

20 SSPE BUFFER. Dissolve 210 g NaCl, 27.6 g NaH₂PO₄, H₂O, 7.4 g Na₂ EDTA, and approximately 8.3 ml of saturated NaOH (to pH 7.0). Adjust volume to one liter, filter, autoclave and store at room temperature.

20 SSC BUFFER. Dissolve 175.3 g NaCl and 88.2 g sodium citrate in 800 ml water. Adjust pH to 7.0 with a few drops of concentrated NaOH. Adjust volume to one liter, filter, autoclave, and store at room temperature.

0.3 M NaOH. Store at room temperature.

1 M and 2 M Ammonium Acetate. Filter and store at room temperature.

4 x BP. Dissolve 2 g each bovine serum albumin, crystalline Pentex Fraction V., and polyvinylpyrrolidone, average mol. wt. = 40,000 in 100 ml water. Filter and store at -20° C in small aliquots.

10% (W/V) Sodium Dodecyl Sulfate (SDS). Store at room temperature.

Other Materials:

Nitrocellulose filters (Schleicher & Schuell, Keene, NH, #BA85) or Zetabind membranes (AMF, Meriden, CT, #NM511-01-04-5SP).

Whatman 3MM chromatography paper.

DEAE-cellulose. Whatman DE52, precycle according to the manufacturer.

Dot blot filtration apparatus (HYBRI•DOT™, BRL, Bethesda, MD; Minifold II, S & S, Keene, NH; BIODOT™, Bio-Rad Laboratories, Richmond, CA.)

Water baths.

Vacuum pump.

Vacuum oven.
Automatic micropipets and tips.
Disposable tuberculin syringes (1 ml).
Glass wool.

5

Sample Preparation

The method of fixing target DNA to nitrocellulose filters will vary depending on the source of the sample. Purified DNA may be bound by filtration by Southern blotting, by Dot blotting. Virus-infected whole
10 cells can be fixed to nitrocellulose by spotting followed by denaturation.

The following procedure is suitable for binding the DNA to nitrocellulose filters. It should be modified according to the needs of the individual experiment. The following steps are performed at room temperature unless indicated otherwise.

NOTE: Nitrocellulose and nylon filters should be stored in a dry cool place and protected from fumes of
15 organic solvents. Handle the filters with gloves.

1. Cut the nitrocellulose or nylon to fit the dot blot apparatus and wet the filters in hot (approximately 60 °C) sterile water. Make certain that the sheet is completely wet, then rinse it several additional times with hot clean water.
2. Soak the sheet in 6 x SSC buffer while preparing the samples.
- 20 3. Resuspend DNA in 0.1 ml 0.3 M NaOH and boil for 5 minutes.
4. Immediately chill the sample in ice and add 0.1 ml 2 M Ammonium Acetate. Mix well.
5. Clamp sheet in the dot blot apparatus, apply a gentle vacuum and wash the wells with 200µl 6 x SSC.
6. Decrease the vacuum and add the sample. A sample should take several minutes to filter in order for the DNA to bind quantitatively.
- 25 7. After the samples have filtered through, wash the wells as in Step 5.
8. Remove the filter from the apparatus, wash it in 6 x SSC, and air dry.
9. Place the filter between the folded sheets of Whatman 3MM paper and bake in a vacuum oven for one hour at 80 °C. Nylon support should be dried for 15-30 minutes at 80 °C.
10. Filter sheets prepared thus may be stored dry at 4 °C in sealed plastic bags.

30

Hybridization Detection

1. Prehybridize the filter in 5 x SSPE, 1 x BP, and 1% SDS for 5 minutes at 50 ° ± 2 °C.
2. Transfer the filter to a boilable polyethylene bag, add 50µl per cm² of 5 x SSPE, 1 x BP, 1% SDS,
35 containing ~ 5-500 ng/ml of appropriately labelled probe and heat-seal the bag. Incubate for 30 minutes at 50 ° ± 2 °C with gentle agitation.
3. After the incubation, open the bag and discard the buffer into an appropriate waste container. Wash the filter as follows, using a large excess of 1 x SSPE, 1% SDS; three times at 37 °C for five minutes each followed by one wash at 50 °C for one minute.
- 40 4. At this point different procedures can be carried out depending upon the nature of the label and detection system (radioisotope, enzyme/color dye system, etc.). With regard to the further teaching in this patent application and in reference to fluorescently labeled probes in particular, the filters would be washed one time with 1 x SSPE. At this point procedures discussed later in the application apply.

The probe is designed and selected to provide a complementary sequence for binding to the target
45 DNA or RNA sequences. Short-chain probes are preferred, such as those containing from 10 to 50 nucleotide units, and a single complementary sequence to a characterizing sequence of the target DNA or RNA. Longer probes such as probes containing in excess of 200 nucleic acid units may contain a plurality of complementary sequences. Such long probes may be used to obtain some of the advantages of the present invention. Probes prepared by nucleotide synthesis which contain from 18 to 30 nucleic acid units
50 are particularly desirable. Methods for synthesizing and labelling such short probes are described in the published PCT Application WO 84/03285.

As in prior practice, any one of a number of different kinds of reporter groups or labels may be attached to the probes. Radioactive labels which may be employed include ³²P, ³H, ¹⁴C labelled probes. However, labels generating some fluorescent or colorimetric response are preferred, as described in the published
55 European Patent Application 0 070 687. In particular, it is preferred to employ fluorescent labels. Fluorescent labelling is described in PCT Publication No. WO 84/03285. The fluorescent labels may include fluorescein, Texas Red, Lucifer Yellow, pyrene, lanthanide complexes, etc. The preferred probe is one containing a single fluorescent label attached to a complementary sequence of 18 to 30 nucleic acid units.

However, the use of two or more fluorophore labels per probe is not excluded, particularly in cases where fluorophore combinations effect a condition for improved detection of fluorescent signal.

The novel steps of the present invention begin following the hybridization and removal of the unbound probe from the initial support matrix (nitrocellulose, etc.). In the event that the concentration of the hybridized probe is sufficient to provide appropriate selectivity and sensitive detection, the method of the present invention need not be employed.

The first step in utilization of the new methodology involves the dehybridization of the labelled probe from the initial support matrix. Dehybridization may be effected by known denaturation procedures, which cause double-stranded nucleic acids to separate into single-stranded form. For example, dehybridization may be carried out by heating the support matrix-target DNA/hybridized labelled probe in a small amount of water or dilute buffer (0.1 to 1X SSPE or SSC) to a temperature sufficient to produce denaturation. In general, effective denaturing temperatures (also called melting temperatures) are in the range from 60° to 90° C. The required temperature depends on the length and "GC" (Guanine & Cytosine) content of the probe. In an alternative procedure, an appropriate buffered solution containing a denaturing agent such as urea or formamide can be used to cause dehybridization of the probe. For example, a 6 molar urea (0.1X - 1X SSPE or SSC) solution is effective in dehybridizing the probe at room temperature (20°-30° C). A variety of general procedures and agents exist for causing dehybridization of the probe from the support. In general, selection is based on the ease and speed of the procedure and the important requirement that it not interfere or inhibit the second stage of the process, i.e., where the labelled probe is adsorbed and concentrated onto a suitable material for detection.

In an alternative procedure which is not known to have been previously employed, the dehybridization of the labelled probe is effected by means of a "displacer probe." The displacer comprises a single-stranded oligonucleotide sequence including the same sequence as the labelled probe with an additional sequence on one side, the additional sequence being complementary to a corresponding sequence on the target strand to which the probe is bound. The additional sequence should be of sufficient length to enable the displacer to preferentially bind to the target sequence near the probe and then efficiently displace the probe. For example, the additional sequences may contain from 10 to 20 nucleic acid units and may be provided on either side of the labelled probe sequence. By way of more specific example, where the hybridized labelled-probe contains from 18 to 30 nucleic acid units, the displacer may correspondingly contain from about 28 to 50 nucleic acid units. Longer displacer probes can also be employed.

In using the displacer, it may be applied in aqueous buffered solution to the hybridized probe. For example, the supported hybridized probe may be immersed in a solution of the displacer. The typical concentrations will range from about 1×10^{-13} to 1×10^{-11} moles of the displacer per 50 microliters of solution (0.1 to 5.0 x SSPE or SSC). The displacement may be carried out at room temperatures (viz. 20-30° C) in less than 30 minutes. In general, the conditions employed are optimized to allow efficient hybridization of the displacer probe/dehybridization of the labelled probe, and as minimal leaching of non-specific bound labelled probe as possible.

This displacement procedure has a further important advantage in that it increases selectivity. During hybridization some of the labelled probe may have become bound to components of the sample other than the target sequence. Where such extraneous binding has occurred, a false positive indication can be obtained if the non-specifically bound labelled probe is detected in the sample. By using the displacer, however, displacement occurs for only that portion of the probe specifically hybridized to the target sequence. Thus, the displacer probe gives double specificity to the assay procedure.

With the preferred procedures described above, the bound labelled probe can be dehybridized from the support material and obtained in the form of an aqueous solution. Direct detection of the labelled probe in the aqueous solution is usually not feasible because the concentration of the labelled probe is extremely low. However, by concentrating and focusing the dehybridized probe, it can be effectively detected at concentrations down to 10^{-18} moles and below.

In one preferred procedure, the dehybridized probe is adsorbed from the solution onto an anion exchange material which provides positively charged groups for binding the probe. Weak basic anion exchange materials containing substituted amino groups may be employed as well as strongly basic anion exchange materials with quaternary amino groups.

It is presently believed, however, that quaternary amino type materials are the most desirable.


Anion exchange materials suitable for use in the method of this invention are available commercially. Such materials can be prepared from substrates comprising cross-linked dextran or glass. Such materials are usually in the form of porous beads or granules of ultra-fine particle sizes (viz. 50-200 micron diameter). Anion exchange materials based on cross-linked dextran are available from Pharmacia (Uppsala, Sweden). These products are in the form of small, generally spherical beads with varying degrees of porosity.

for ultrasensitive detection. Usually the analysis and quantitation takes only a matter of minutes for a given sample.

The method of this invention can be further understood in relation to the following examples.

5 EXAMPLE I

The following example demonstrates the advantage for detecting fluorescent probes concentrated onto suitable materials, like Sephadex QAE A-25 beads, over observing the same probe bound to nitrocellulose filter material. The fluorescent probe used in these experiments was a 22-mer sequence specific for a
10 genomic target sequence present in Herpes Simplex Virus (HSV). The probe was labeled with a Texas Red fluorophore group having an excitation and emission maxima at 596nm and 620nm respectively. The labeled probe was designated TR-HSV, and the sequence and position of the label are shown below:

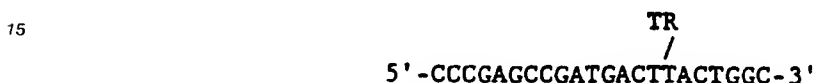

The first experiment involved the epifluorescent analysis of the TR-HSV probe bound to nitrocellulose
20 filters. A 10-fold dilution series of TR-HSV probe ranging from 2×10^{-12} to 2×10^{-18} mole per microliter of 1X SSPE buffer was prepared. One microliter of each dilution was adsorbed onto nitrocellulose filter material forming a 2-3 millimeter diameter spot. A control spot of 1 μ L of 1 x SSPE (no probe) was also prepared. The spots were now analyzed by using a photon counting epifluorescent microscope filtered for Texas Red fluorescence (Ex ~ 550nm; Em ~ 630nm). The fluorescent signal from each sample spot was
25 counted (EG&G Photon Counting System) for one second (each sample spot was counted twice). Results for the analysis are shown in Table I. Since only 1/10 the area of each spot was analyzed (10x microscope field), samples 1 through 7 range from 2×10^{-13} to 2×10^{-19} mole in actual TR-HSV probe concentration.

30

TABLE I
TR-HSV Probe/Nitrocellulose Filter

35	<u>Sample</u>	<u>TR-HSV (mole)*</u>	<u>Avg. Counts (1 second)</u>
	1	2×10^{-13}	1.5×10^6
	2	2×10^{-14}	157,429
40	3	2×10^{-15}	27,918
	4	2×10^{-16}	26,246
	5	2×10^{-17}	21,155
45	6	2×10^{-18}	22,100
	7	2×10^{-19}	23,397
	Control	---	22,399

50

*Concentration in area observed

The results in Table I show that fluorescence from the TR-HSV probe can be detected down to about 2×10^{-16} mole (sample 4) on nitrocellulose filter material.
55

The second experiment involved the epifluorescent analysis of the TR-HSV probe bound to QAE-Sephadex A-25 beads. A 10-fold dilution series of TR-HSV probe ranging from 4×10^{-11} to 4×10^{-17} mole per 10 microliters of 1 X SSPE buffer was prepared. Approximately 1000 QAE-Sephadex A-25 beads

(previously swollen in 1 X SSPE buffer) were added to each of the sample dilutions. A control sample of beads with no TR-HSV probe was also prepared. The samples were allowed to stand for about one hour with intermittent gentle agitation. Samples were then analyzed using the photon counting epifluorescent microscope. In this experiment, the fluorescence from individual beads in each sample were counted. The results in Table 2 are the average counts (one second) from two beads in each sample counted twice. The results are given in terms of TR-HSV concentration per bead, i.e., it is assumed most of the probe has adsorbed evenly to all the beads in the sample.

TABLE 2

TR-HSV Probe/QAE-Sephadex A-25

Sample	TR-HSV (mole)/bead	Avg. Counts (1 second)
1	4×10^{-14}	9×10^6
2	4×10^{-15}	645,867
3	4×10^{-16}	43,229
4	4×10^{-17}	15,614
5	4×10^{-18}	13,581
6	4×10^{-19}	5,429
7	4×10^{-20}	4,782
Control	---	4,320

The results in Table 2 show that fluorescence from the TR-HSV probe at 4×10^{-20} mole per bead still has counts significantly higher than the control beads (no probe).

The third experiment is basically a complete repeat of the second experiment. The results are given in Table 3.

TABLE 3

TR-HSV Probe/QAE Sephadex A-25

Sample	TR-HSV (mole)/bead	Avg. Counts (1 second)
1	4×10^{-14}	1×10^7
2	4×10^{-15}	1×10^6
3	4×10^{-16}	85,939
4	4×10^{-17}	17,022
5	4×10^{-18}	13,627
6	4×10^{-19}	10,182
7	4×10^{-20}	8,082
Control	---	7,341

The results in Table 3 show again the fluorescence from the TR-HSV probe at 4×10^{-20} mole per bead has counts significantly (~10%) higher than the control beads. The higher counts for control beads in experiment 3 (7,341) compared to experiment 2 (4,320) are due to the higher excitation intensity used in

this experiment.

In summary, TR-HSV probe can be detected down to about 4×10^{-20} mole on a single bead of QAE-Sephadex A-25 by epifluorescence analysis. On nitrocellulose filter material it can be detected down to only about 2×10^{-16} mole.

EXAMPLE II

In further experiments similar to those described in Example 1 it was demonstrated that another type of fluorescent probe (Fluorescein labelled) can be concentrated and more favorably detected on QAE-Sephadex A-25 beads. Also, it is shown that other supports, such as QAE Controlled Pore Glass Beads, work favorably.

The first experiment involves the epifluorescent analysis of a fluorescein labelled HSV probe (F-HSV) bound to Zeta Bind (nylon) filter material. The F-HSV probe is exactly the same as the TR-HSV probe except for the different fluorophore label. The fluorescein label has its excitation and emission maximum at ~490nm and ~520nm respectively. The experiment with F-HSV probe adsorbed to Zeta Bind was carried out exactly like the first experiment in Example I. However, the Photon counting epifluorescence microscope was filtered for fluorescein fluorescence (Ex ~ 450 nm, Em ~ 520 nm). Results for the analysis are shown in Table 4.

TABLE 4

F-HSV Probe/Zeta Bind Filter

<u>Sample</u>	<u>F-HSV (mole)*</u>	<u>Avg. Counts (1 second)</u>
1	6×10^{-13}	843,568
2	6×10^{-14}	200,248
3	6×10^{-15}	55,914
4	6×10^{-16}	45,587
Control	---	47,675

*Concentration in area observed.

The results from Table 4 show that fluorescence at $\sim 6 \times 10^{-15}$ mole of probe is just distinguishable from the control sample.

The second experiment involves the epifluorescent analysis of the F-HSV probe bound to QAE-Sephadex A-25 beads. This experiment was carried out in a manner similar to the second experiment in Example I. Results for the analysis are shown in Table 5.

TABLE 5

F-HSV Probe/QAE-Sephadex Beads

Sample	F-HSV (mole)/bead	Avg. Counts (1 second)
1	7×10^{-14}	10^7
2	7×10^{-15}	2×10^6
3	7×10^{-16}	136,736
4	7×10^{-17}	56,924
5	7×10^{-18}	31,401
6	7×10^{-19}	32,194
Control	---	27,507

The results from Table 5 shown that fluorescence from the F-HSV probe is detectable at the 7×10^{-18} mole per bead level.

As was the case with the TR-HSV probe (Example I), the F-HSV probe is significantly more detectable on QAE-Sephadex A-25 beads than on Zeta Bind filter material.

The third experiment involves the epifluorescent analysis of TR-HSV probe bound to QAE-Controlled Pore Glass beads. This experiment was carried out similar to the others described above. Results for the analysis are shown in Table 6.

TABLE 6

TR-HSV Probe/QAE-Controlled Pore Glass Beads

Sample	TR-HSV Probe/Bead	Avg. Counts (1 second)
1	4×10^{-14}	10^7
2	4×10^{-15}	2×10^6
3	4×10^{-16}	179,131
4	4×10^{-17}	51,797
5	4×10^{-18}	20,454
6	4×10^{-19}	16,897
7	4×10^{-20}	15,497
Control	---	11,930

The results from Table 6 show that TR-HSV probe is detectable to the 4×10^{-20} mole level on individual QAE-Controlled Pore Glass beads.

EXAMPLE III

The experiment described below demonstrates that fluorescent probe actually hybridized to a target DNA sequence immobilized on a polystyrene bead support material (not suitable for fluorescent analysis) can be dehybridized and concentrated onto QAE-Sephadex A-25 beads for epifluorescent analysis. In this experiment the TR-HSV probe used was also labelled with 32 Phosphorous radioisotope (32 P). Thus a comparison could be made between the fluorescent and radioisotope level of detection.

Using one milligram samples of polystyrene beads (~1-2 micron diameter, substituted with a complementary target sequence), a 10-fold dilution series of ^{32}P -TR-HSV probe ranging from $\sim 1 \times 10^{-13}$ to $\sim 1 \times 10^{-18}$ mole in 50 μl of 5X SSPE hybridization buffer were prepared (control sample contained no probe). The samples were hybridized for 30 minutes at 40°C, then centrifuged and supernatant removed. The samples were then washed three times with 5X SSPE. After the final wash and removal of supernatant, 10 μl of 6M Urea (1X SSPE) was added to each of the samples in order to dehybridize the probe from the support. Dehybridization was carried out for about 30 minutes at room temperature. At this point the polystyrene beads were separated from the solution containing dehybridized ^{32}P -TR-HSV probe by centrifugation. To this solution between 20 to 30 QAE-Sephadex A-25 beads were added. The dehybridized probe was allowed to adsorb to the QAE-Sephadex beads for about 30 minutes. At this point the QAE-Sephadex A-25 beads were separated from the solution. The beads (all 20 to 30) were first counted on a scintillation counter (one minute counts) to determine level of probe finally transferred to the QAE-Sephadex A-25 beads by radioisotope method. After scintillation counting individual beads were counted on the photon counting epifluorescent microscope to determine sensitivity by fluorescence. The results for each type of analysis are shown in Table 7.

TABLE 7
 ^{32}P -Tr-HSV Transfer to QAE-Sephadex

<u>Sample</u>	<u>Total Conc. ^{32}P-Tr-HSV Probe</u>	<u>^{32}P Counts</u>	<u>^{32}P-Tr-HSV Conc. per bead</u>	<u>Fluorescent Counts (1 sec.)</u>
1	$\sim 6 \times 10^{-15}$	7460	$\sim 3 \times 10^{-16}$	104,572
2	$\sim 6 \times 10^{-16}$	951	$\sim 3 \times 10^{-17}$	28,847
3	$\sim 6 \times 10^{-17}$	159	$\sim 3 \times 10^{-18}$	13,704
4	$\sim 6 \times 10^{-18*}$	43	$\sim 3 \times 10^{-19}$	13,170
5	$\sim 6 \times 10^{-19*}$	40	$\sim 3 \times 10^{-20}$	10,227
Control	---	44	---	8,900

*extrapolated

The results in Table 7 show that about 6×10^{-17} mole of ^{32}P -TR-HSV probe transferred to QAE-Sephadex beads could be detected by radioisotope counting (^{32}P , 1 minute count) before control or background level is reached (~44 counts). Fluorescent counts were found to be above the control through the whole sample series 1→5. When considered in terms of individual beads, on the order of 10^{-20} mole of ^{32}P -TR-HSV probe was detectable.

EXAMPLE IV

The following experiment demonstrates the usefulness of a so-called "displacer probe" in effecting the dehybridization of a fluorescence labelled probe hybridized to target DNA immobilized on a support material.

In this experiment ^{32}P -TR-HSV probe was hybridized to plasmid target DNA (pHSV106) immobilized onto nitrocellulose filter material. The pHSV106 plasmid DNA was immobilized onto the nitrocellulose filter using the Hybrid-Dot apparatus and experimental procedure described earlier in this patent application. The concentrations of pHSV106 plasmid DNA in each "dot" (~2 millimeter diameter spot) were 1) 100ng; 2) 10 ng, 3) 1 ng, 4) 100 pg, and 5) control (non-target DNA). ^{32}P -TR-HSV probe was hybridized to the pHSV106 target dots by the procedures discussed earlier in this application. After hybridization and appropriate washing, the level of ^{32}P -TR-HSV probe hybridized to each "dot" was measured by a scintillation counter. The results are given in Table 8. After counting, the dehybridization procedure was carried out using a specific displacer probe for the ^{32}P -TR-HSV probe. The sequence of the displacer probe is given below:

5'-CCCCCAGCACCTGCCAGTAAGTCATCGGCTCGGG-3'

The portion of the displacer probe underlined is identical to that of the ³²P-TR-HSV probe. Dehybridization of ³²P-TR-HSV/pHSV106 target DNA nitrocellulose "dots" was carried out by treating each one of the sample "dots" with 50μl of 1X SSC buffer solution containing approximately 100ng of the displacer probe for 30 minutes at ~25° C. The nitrocellulose filter "dots" were then separated from the solution and each solution and "dot" sample was counted on the scintillation counter. The results are given in Table 8.

TABLE 8

Displacer Probe Experiment

Sample (pHSV106)	1st Count (1 min.) after Hybridization & Wash	2nd counts (1 min) After Displacer Probe	
		Filter "Dots"	Solution
1 100ng	39,672	5,969	32,353
2 10ng	5,805	1,347	4,296
3 1ng	1,702	1,113	754
4 100pg	1,543	1,109	706
5 None	1,296	1,210	198

The results in Table 8 show that the "displacer probe" treatment has successfully dehybridized the ³²P-TR-HSV probe from the pHSV106 target/nitrocellulose filter "dots."

Claims

1. The hybridization assay method for detecting a specific polynucleotide target sequence wherein a nucleic acid containing target sample in single stranded form is affixed to a support, contacted thereon under hybridizing conditions with a reporter group-labeled single-stranded polynucleotide probe having a sequence complementary to the target sequence, at least part of said probe being bound to said affixed sample when the target is present, said probe hybridizing with said target sequence, then removing the unbound portion of said probe from said support, and detecting the presence of said probe by means of its reporter group, wherein the sensitivity of the assay is increased by employing the additional steps comprising:
 - (a) after removal of the unbound portion of the probe, dehybridizing the bound portion;
 - (b) separating the dehybridized portion from the support and the residual sample; and
 - (c) concentrating the dehybridized separated probe before detecting its presence by means of said reporter group.
2. The method of claim 1 in which the dehybridized portion of the probe is concentrated by adsorption from an aqueous solution onto a basic anion exchange material binding said probe.
3. The method of claim 2 in which the anion exchanger functional groups are quaternary amine groups.
4. The method of claim 3 in which said anion exchange functional groups are provided on a substrate selected from cross-linked dextran or glass.
5. The method of claim 1 in which the bound labelled probe portion is dehybridized by contacting the supported sample with a displacer probe comprising a single stranded polynucleotide sequence including the same sequence as the labelled probe with an additional sequence on at least one side thereof complementary to the corresponding sequence on the target sequence to which the labelled probe is bound, said additional sequence being of sufficient length to cause said displacer to preferentially bind to said target sequence and subsequently displace said labelled probe therefrom.

6. The method of claim 1 in which said probe comprises a sequence of 10 to 50 nucleotide units and has a reporter group attached thereto, and said probe is concentrated on said anion exchange material before detecting its presence by means of said reporter groups.
- 5 7. The method of claim 6 in which said anion exchange material comprises an adsorbent substrate selected from cross-linked dextran or glass with substituted amino groups attached thereto, and is in the form of generally spherical beads.
8. The method of claim 6 in which said probe is labelled with a fluorescent reporter group.
- 10 9. The method of claim 6 in which the bound probe portion is dehybridized by contacting the supported sample with a displacer comprising a single-stranded polynucleotide sequence including the same sequence as the probe with an additional sequence on at least one side thereof complementary to the corresponding sequence on the same side of the target sequence to which the probe is bound, said
15 additional sequence being of sufficient length to cause said displacer to preferentially bind to said target sequence and displace said probe therefrom.

Revendications

- 20 1. Procédé d'essai d'hybridation pour détecter une séquence polynucléotidique-cible spécifique, selon lequel on fixe sur un support, un échantillon-cible contenant des acides nucléiques sous forme monocaténaire, on le met en contact, dans des conditions d'hybridation, avec une sonde polynucléotidique monocaténaire marquée avec un groupe indicateur, ayant une séquence complémentaire de la
25 séquence-cible, au moins une fraction de la sonde étant liée à l'échantillon fixé lorsque la cible est présente, la sonde s'hybride avec la séquence-cible, on élimine ensuite du support la fraction non liée de la sonde, et on détecte la présence de la sonde par l'intermédiaire de son groupe indicateur, la sensibilité de l'essai étant accrue en effectuant les opérations supplémentaires consistant:
 - (a) à déshybrider la fraction liée, après avoir éliminé la fraction non liée de la sonde;
 - (b) à séparer la fraction déshybridée du support et de l'échantillon résiduel; et
 - 30 (c) à concentrer la sonde déshybridée séparée, avant de détecter sa présence par l'intermédiaire du groupe indicateur.
2. Procédé selon la revendication 1, dans lequel la fraction déshybridée de la sonde est concentrée par adsorption à partir d'une solution aqueuse, sur un agent d'échange d'anion basique, se liant avec la
35 sonde.
3. Procédé selon la revendication 2, dans lequel les groupes fonctionnels échangeurs d'anion sont des groupes amine quaternaires.
- 40 4. Procédé selon la revendication 3, dans lequel les groupes fonctionnels échangeurs d'anion, sont sur un substrat choisi parmi le dextrane réticulé et le verre.
5. Procédé selon la revendication 1, dans lequel la fraction liée de la sonde marquée est déshybridée en mettant en contact l'échantillon associé au support, avec une sonde de déplacement comprenant une
45 séquence polynucléotidique monocaténaire contenant la même séquence que la sonde marquée et une séquence supplémentaire à au moins l'une de ses extrémités, complémentaire de la séquence correspondante de la séquence-cible à laquelle la sonde marquée est liée, la séquence supplémentaire ayant une longueur suffisante pour que ladite sonde de déplacement se lie de préférence à la séquence-cible, et se sépare ensuite la sonde marquée de celle-ci.
- 50 6. Procédé selon la revendication 1, dans lequel la sonde comprend une séquence de 10 à 50 motifs nucléotidiques, et comporte un groupe indicateur lié, et on concentre la sonde sur l'agent d'échange d'anion, avant de détecter sa présence par l'intermédiaire des groupes indicateurs.
- 55 7. Procédé selon la revendication 6, dans lequel l'agent d'échange d'anion comprend un substrat adsorbant choisi parmi le dextrane réticulé ou le verre auquel des groupes amino substitués sont liés, et est sous la forme de perles globalement sphériques.

8. Procédé selon la revendication 6, dans lequel la sonde est marquée avec un groupe indicateur fluorescent.
- 5 9. Procédé selon la revendication 6, dans lequel la fraction de sonde liée est déshybridisée, en mettant en contact l'échantillon associé au support, avec un agent de déplacement comprenant une séquence polynucléotidique monocaténaire contenant la même séquence que la sonde, avec une séquence supplémentaire à au moins l'une de ses extrémités, complémentaire de la séquence correspondante à la même extrémité de la séquence-cible à laquelle la sonde est liée, la séquence supplémentaire ayant une longueur suffisante pour que ledit agent de déplacement se lie de préférence à la séquence-cible, et sépare la sonde marquée de celle-ci.
- 10

Patentansprüche

- 15 1. Das Hybridisierungstestverfahren zum Nachweisen einer spezifischen Polynukleotidzielsequenz, worin eine Nukleinsäure enthaltende Zielanalysenprobe in einzelsträngiger Form auf einem Träger fixiert ist, die darauf unter Hybridisierungsbedingungen mit einer mit einer Reportergruppe markierten einzelsträngigen Polynukleotidprobe in Kontakt gebracht wird, die eine zu der Zielsequenz komplementäre Sequenz besitzt, wobei zumindest ein Teil der Probe an die fixierte Analysenprobe gebunden wird, wenn das Ziel anwesend ist, wobei die Probe mit der Zielsequenz hybridisiert, dann Entfernen des ungebundenen Anteils der Probe von dem Träger, und Nachweisen der Anwesenheit der Probe anhand seiner Reportergruppe, worin die Empfindlichkeit des Tests gesteigert wird durch Anwendung der zusätzlichen Schritte umfassend:
 - a) nach Entfernen des ungebundenen Anteils der Probe, Enthybridisieren des gebundenen Anteils;
 - 25 b) Abtrennen des enthybridisierten Anteils von dem Träger und der restliche Analysenprobe; und
 - c) Konzentrieren der enthybridisierten abgetrennten Probe vor dem Nachweis ihrer Anwesenheit anhand der Reportergruppe.
2. Das Verfahren nach Anspruch 1, in dem der enthybridisierte Anteil der Probe konzentriert wird durch Adsorption aus einer wässrigen Lösung an ein basisches Anionenaustauschermaterial, das die Probe bindet.
- 30 3. Das Verfahren nach Anspruch 2, in dem die funktionellen Gruppen des Anionenaustauschers quartäre Aminogruppen sind.
- 35 4. Das Verfahren nach Anspruch 3, in dem die funktionellen Gruppen des Anionenaustauschers auf einem Substrat bereitgestellt werden, ausgewählt aus vernetztem Dextran oder Glas.
5. Das Verfahren nach Anspruch 1, in dem der gebundene, markierte Probenanteil enthybridisiert wird durch In-Kontakt-Bringen der Probe auf dem Träger mit einer verdrängenden Probe, umfassend eine 40 einzelsträngige Polynukleotidsequenz, die die gleiche Sequenz wie die markierte Probe enthält mit einer zusätzlichen Sequenz an mindestens einer Seite davon, die komplementär ist zu der entsprechenden Sequenz der Zielsequenz, an die die markierte Probe gebunden ist, wobei die zusätzliche Sequenz von ausreichender Länge ist, um den Verdränger zu veranlassen, bevorzugt an die Zielsequenz zu binden und anschließend die markierte Probe davon zu verdrängen.
- 45 6. Das Verfahren nach Anspruch 1, in dem die Probe eine Sequenz von 10 bis 50 Nukleotideinheiten umfaßt und eine daran angefügte Reportergruppe besitzt, und die Probe auf dem Anionenaustauschermaterial konzentriert wird, bevor ihre Anwesenheit anhand der Reportergruppen nachgewiesen wird.
- 50 7. Das Verfahren nach Anspruch 6, in dem das Anionenaustauschermaterial ein Adsorptionssubstrat umfaßt, ausgewählt aus vernetztem Dextran oder Glas mit daran angefügten substituierten Aminogruppen, und es in der Form von normalerweise runden Kügelchen vorliegt.
8. Das Verfahren nach Anspruch 6, in dem die Probe mit einer fluoreszierenden Reportergruppe markiert ist.
- 55 9. Das Verfahren nach Anspruch 6, in dem der gebundene Probenanteil enthybridisiert wird durch In-Kontakt-Bringen der Probe auf dem Träger mit einem Verdränger, umfassend eine einzelsträngige

Polynukleotidsequenz, die die gleiche Sequenz wie die Probe enthält mit einer zusätzlichen Sequenz an mindestens einer Seite davon, die komplementär zu der entsprechenden Sequenz auf der gleichen Seite der Zielsequenz ist, an die die Probe gebunden ist, wobei die zusätzliche Sequenz von ausreichender Länge ist, um den Verdränger zu veranlassen, bevorzugt an die Zielsequenz zu binden und die Probe davon zu verdrängen.

5

10

15

20

25

30

35

40

45

50

55